

permobil

TILITE MANUAL WHEELCHAIRS

2GX

TITANIUM PERFORMANCE

The 2GX combines the advantages of a modular folding frame with the high strength and vibration dampening characteristics of titanium. Add full adjustability and the 2GX is designed to empower - no matter where you choose to roll.





TECHNICAL SPECIFICATIONS

Frame Features	Titanium Box Frame Folding Modular Design Full Adjustability	Camber	0°, 2°, 4°, or 6°
Total Width Folded*	10.75" (27.3cm)	Centre of Gravity Adjustment	6" (15.2cm)
Total Chair Width*	23" (58.4cm)	Transit Options	Transport Bracket
Seat Width	12" – 22" (30.4 – 55.8cm)	Transport Weight*	Fixed Front - 5.9 kg. Swing Away - 7 kg.
Seat Depth	12" – 20" (30.4 – 50.8cm)	Weight Limit	120.2 kg Standard 158.7 kg Heavy Duty
Front Seat Height	13.5" – 21" (34.2 – 53.3cm)	Warranty	Five Years on Frame One Year on TiLite Components
Rear Seat Height	13" – 20" (33 – 50.8cm)		
Front Frame Angle	Fixed Front - 70°, 80°, or 85° Swing Away - 70°, 80°, 85°, or 90°		
Seat Back Type	<i>Folding Seat Back</i> Adj. Height 11" – 20.5" (27.9 - 52cm) Adj. Angle 80° – 101° <i>Fixed Seat Back</i> Adj. Height 11" – 20" (27.9 - 50.8cm)		

* Weights and widths based on 16" x 16" chair with 0° camber, minimum rear wheel spacing and all standard features. Rear wheels are not included in transport weight. For reference, a pair of standard 24" rear wheels weighs 3.2 kg.

MODULAR DESIGN


Optimise your ride with custom configuration. This chair's modular design allows you to add options including after market seating systems, switching fixed and swing away front ends, and growing chair width and / or depth to meet changing customer needs.

CUSTOM CONFIGURATION

Custom Configuration allows you to add options and components to your chair to meet your needs and to optimise your riding experience. This includes choices of wheels, handrims, tyres, wheel locks, colour, tattoos, and more. TiLite offers more options than any other manufacturer!


COLOUR OPTIONS

HIGH GLOSS PAINT

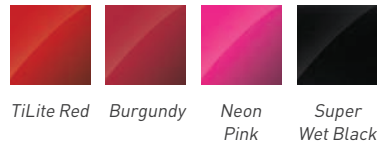

TiLite Red Burgundy Neon Pink Super Wet Black

MATTE PAINT

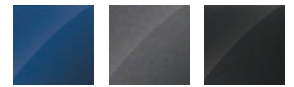

Deep Blue Ice Titanium Grey Black Lava

METALLIC PAINT



Salsa Red Metallic Tangerine Metallic Acid Green Metallic Ocean Blue Metallic

CANDY PAINT



Ruby Red Azure Blue Electric Plum

NATURAL



Satin Titanium Polished Titanium

PEARLESCENT PAINT



Midnight Blue Pearl White Pearl