

permobil

Bariatric Collection



Bed Collection

Botero Bariatric Bed

The specific bariatric design focus offers width and length adjustable deck to manage changing needs in both ward and respite environments.

- Low height of 38cm.
- Antibacterial powder coating for improved infection control management.
- Backrest translation (autoregression) reduces compression of the thoracic and abdominal areas as the backrest is raised and alleviates shear and pressure for improved skin integrity
- Compatible with all leading bed movers.
- Attendant control panel with intuitive pre-programmed single touch functions for speed and simplicity.
- 15cm TENTE castors.
- Central locking castors with integrated directional lock.
- Full Trendelenburg functionality.
- Integrated folding side rails.
- Removable width adjustable head and foot panels.
- Back rest angle 0 to 70 degrees.
- Knee bend angle 0 to 30 degrees.
- Trendelenburg -17 to 17 degrees.

SWL 400Kg

Specifications

Size	Dimensions WxLxH
Width	102cm - 122cm
Length	226cm - 244cm
Height	38cm - 85cm



Deutscher Bariatric Floorline Bed - Big Ted

The all new Deutscher Bariatric Bed has a safe working load rating of 350Kg and has been released onto the market after many months of stringent testing to ensure it meets the requirements necessary for the safe handling of larger clients.

The Big Ted Bed is a 3 function bed with a generous mattress platform of 107cm x 203cm.

It is fitted with twin actuators on both the back rest and the knee bend functions and a pair of strong actuators at each end of the bed to raise and lower the mattress platform. A larger capacity controller and transformer is fitted to handle the power. This new floorline bed uses Dewert Electrics, which are well represented and supported in Australia, meaning spare parts are readily available, should they be required.

The height range of the bed is 13.5cm to 73.5cm.

The bed comes standard with a hand controller as well as a built in control panel in the foot end tower.

A pump bracket is fitted to the foot end tower to accommodate an alternating mattress, if required.

For ease of mobility, heavy duty dual 12.5cm castor wheels have been fitted with quick action locking bars at each end of the bed.

The ability to be easily folded when not in use is an enormous advantage for moving in tight areas, or for when placing the unit into storage.

Available in King Single

SWL 350Kg

Specifications

Size	Dimensions WxLxH
Width	107cm
Length	203cm
Height	135cm - 760cm



Mattress Collection

Aspire Lifecomfort Active Air BARI 10 Mattress Replacement

Integrated Turning Therapy. Allows patient rotation up to 400 for improved pressure care and to assist staff with patient handling requirements.

Rapid Pump Inflation. Maximises therapy response times and minimises staff resources by fully inflating in only 2x minutes.

Hospital Grade Cover. Thicker weight polyurethane fabric with increased bleach resistance to handle modern hospital infection control standards. Suitable for Sodium Hydrochloride solutions.

Adaptable Pump. Pump can be switched between modes for typical users (up to 160Kg) and bariatric users (up to 450Kg) ensuring a single asset can meet the broadest spectrum of patient needs.

SWL 30Kg - 450Kg

Specifications

Size	Dimensions WxLxH
King Single	107cm x 200cm x 27cm
Super King Single	123cm x 200cm x 27cm



Talley QUATTRO® Sentinel Mattress Replacement

Suitable for the prevention and treatment of pressure ulcers within all care environments and is specifically designed for the care of bariatric patients.

The QUATTRO Sentinel mattress is able to deliver safe, effective pressure area care for the prevention and management of pressure related tissue damage in high risk bariatric patients.

100% of staff reported that the QUATTRO Sentinel mattress was effective and very effective at redistributing pressure and maintaining patients skin integrity and was easy to very easy to use.

SWL 370Kg

Specifications

Size	120cm (86 cm without bolsters) x 200cm
Height	18cm
Weight Limit	370Kg
Number of Cells	17
Alternating Cycle	16 minutes



Active 1-in-4 Cell Cycle

370Kg Weight Limit

Removable Side Bolster

Ortho-Differential Support™

Talley Quattro Acute Mattress Replacement

Suitable for the prevention and treatment of pressure ulcers within all care environments, the QUATTRO ACUTE mattress replacement system features Active QUATTRO THERAPY and continuous low pressure therapy modes, together with AUTO support fully automatic pressure adjustment and the unprecedented low interface pressures of DEEP CELL THERAPY.

SWL 250Kg

Specifications

Size	198cm x 1100cm
Height	24cm
Weight Limit	250Kg
Number of Cells	27
Alternating Cycle	Variable



Active 1-in-4 Cell Cycle

Tissuegard™

Deep Cell Therapy™

Ortho-Differential Support™

Modular Therapy Mattress (MTM)

Excellent Pressure Redistribution is achieved by using a combination of foam and ROHO®.

The use of ROHO Overlay inserts minimise risk from internal shearing and external friction. Transfers are made easily as mattress is only 12.5cm in height, same as a standard mattress.

Mattress inner is fully protected by a two way stretch, vapour permeable, incontinent cover with welded seams, this fully fitted mattress cover complies with Hospital Infection Control protocols.

No external power sources required means NO pumps running, soothing environment, NO agitating "white noise" which allows patients to sleep.

Immersion surface gives stable base support for patients unlike unstable alternating surfaces.


Modular system allows high level resources to be utilised only where required and lowers the costs of overall equipment pool.

Mattress is designed specifically to reduce shearing and friction issues created with electric beds allowing clients to sit semi-elevated for long periods.

SWL No Weight Limit

Specifications

Size	Dimensions WxLxH
King Single	107cm x 202cm x 12.5cm

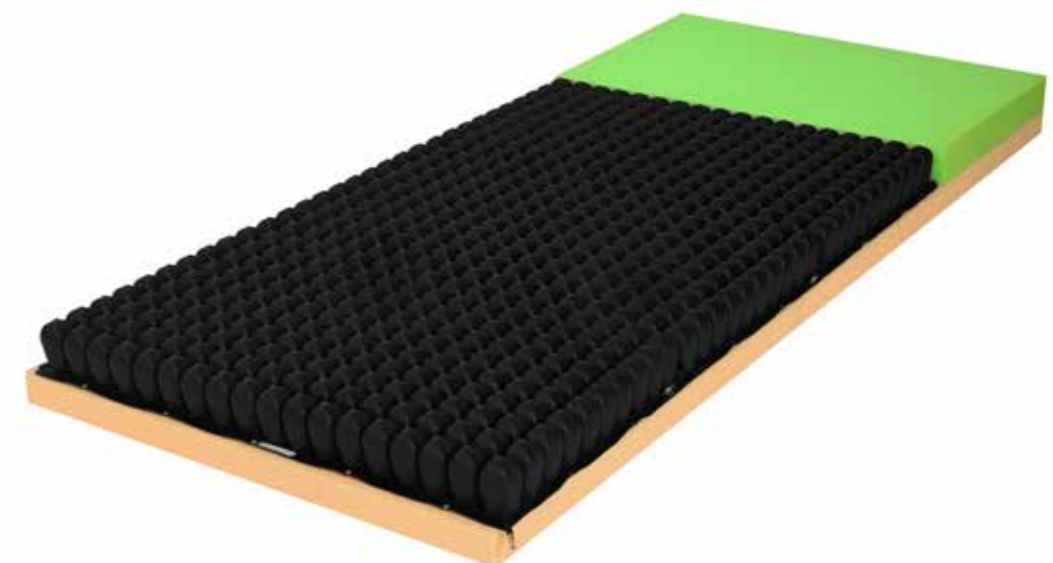

Bathroom and Commode Collection

Aspire MAXI Shower Stool

Super heavy duty adjustable width shower stool. Soft padded seat and arms provide postural support.

Reinforced aluminium frame accommodates bariatric users and is corrosion resistant.

Height adjustable legs with easy clip adjustment points. Raised armrests to assist the user in transferring on and off the stool.

Durable non slip feet provide excellent grip and stability in wet conditions. Stainless steel washer enhances durability.

Size Width 57.5cm - 77.5cm x Depth 45.5cm
SWL 310Kg



Aspire MAXI Over Toilet Aid

Super heavy duty adjustable width over toilet aid. Soft padded flip up seat for easy cleaning.

Reinforced aluminium frame accommodates bariatric users and is corrosion resistant.

Height adjustable legs with easy clip adjustment points. Raised armrests to assist the user in transferring on and off the stool.

Durable non slip feet provide excellent grip and stability in wet conditions. Stainless steel washer enhances durability.

Size Width 57.5cm - 77.5cm x Depth 45.5cm
SWL 310Kg



K Care Shower Commode Maxi

The chair has simple to use swing up arms with padded arm rests and adjustable height swing away foot rests, a back rest with 3 flexible straps for comfort, an integrated push handle and 12.5cm braking castors, make the unit easy to manoeuvre and secure when transferring.

The shower commode is durable and long lasting with a stainless steel frame, impervious PVC straps and water resistant castors.

Size Width 65cm x Depth 46cm
SWL 200Kg



Aspire 600mm Weight Bearing Platform Shower Commode

Heavy duty reinforced wide frame can support bariatric users up to 400kg.

4 x total locking stainless steel castors provide brake access from all angles.

Height adjustable reinforced armrests offer drop flat positioning to assist with lateral slide transfers.

Easy slide weight bearing platform footrest slides completely away for transfers.

Size Width 60cm x Depth 48cm
SWL 400Kg



Aspire Bariatric Shower Commode

Super heavy duty bariatric design with 79cm seat. Anti-bacterial frame coating for improved infection control.

Removable arms to assist with positioning and allow lateral slide transfers.

Height adjustable, swing away footrests for optimal positioning and ease of transfers.

Padded seat and arms for maximum comfort.

Size Width 79cm x Depth 48cm
SWL 350Kg



Aspire Tilt and Recline Shower Commode

30° of backrest angle adjustment in 10° increments for effortless repositioning and maximum user comfort.

-5° to 35° of gas assisted tilt-in-space functionality minimises manual handling efforts during the showering process.

4 x total locking stainless steel castors provide brake access from all angles.

Size Width 46cm - 53cm x Depth 53cm - 58cm
SWL 250Kg



Chair Collection

Aspire Maxi Day Chair

Seat with anterior tilt. Facilitates standing from sitting and provides greater comfort for stiff or painful hips.

Height Adjustable. Easy push button adjustment for comfortable positioning and changing needs.

Depth and backrest height adjustable with ease.

Push handle and transport wheels. Allows the carer to tilt the chair backwards on the rear leg wheels, for ease of transport (without a patient in the chair).

Size Width 60cm - 81cm x Depth 39cm - 53cm
SWL 310Kg



Greiner Bariatric Chair Multiline XXL

Uniquely developed for bariatric patients. An ideal solution in the bariatric space for patient transport, examination and treatment, post-op recovery, bariatric surgery.

Reliable and sturdy mechanic construction with comfortable upholstery.

All adjustment is electronically controlled for maximum patient comfort and independence, as well as reduced nursing burden.

Capable of multiple positioning, including a lay flat position.

Size Width 70cm x Depth 57cm (Overall Depth 110cm to 130cm)
SWL 300Kg



K Care Shower Commode Maxi Tilt (Electric)

Assists transportation and positioning of the larger user for showering and can be positioned over the toilet for toileting.

Electric tilt-in-space mechanism allows carer to adjust the vertical position from 0 to 30 (degrees) ensuring user comfortably retains correct seating position.

Large calf pads to support the legs.

Large backrest with headrest extension to support the upper body and head.

Size Width 60cm x Depth 47cm
SWL 225Kg



K Care Maxi Deluxe V2 Cloud Chair

The K Care Maxi Deluxe V2 Cloud Chair is the market leader in pressure care seating, with unparalleled pressure redistribution and adjustable positioning.

The Maxi Deluxe V2 Cloud Chair offers ease of handling for carers and maximum comfort for patients.

Unique Air Comfort Seating System for total surface pressure reduction for larger users.

Unique power assist tilt-in-space of seat, backrest, legrest and full recline.

Adjustable wings and arms for easy transfer.

Size Width 75cm x Depth 55cm
SWL 300Kg



Cushion Collection

ROHO® Select Series


The QUADTRO SELECT® cushion features four independent compartments of air cells for an individualised fit, each controlled by our exclusive ISOFLO® Memory Control.

ROHO's ISOFLO® Memory Control, located at the front of all Select Series cushions, allows air to flow through the cushion in two ways.

Cushion Collection available in various sizes.

Size Width 20" - 28" x Depth 20" - 28"

ROHO® QUADTRO SELECT®

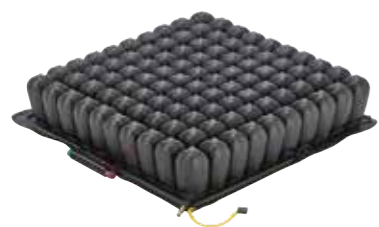

HIGH PROFILE®

The tall interconnected air cells envelop your body, delivering maximum immersion for the ultimate in skin protection, positioning, stability and comfort.



LOW PROFILE®

The interconnected air cells deliver sitting stability. It could be used as a secondary cushion.



ISOFLO Memory Control shown in the OPEN position.

Air can flow freely throughout the entire cushion, similar to a ROHO Single Compartment Cushion.



Power Wheelchair

M300 Corpus® HD

The M300 Corpus® HD, is a heavy duty power wheelchair designed specifically for users up to 204Kg.

The M300 Corpus® HD offers

- A highly maneuverable midwheel drive power base.
- Power tilt (and power recline and power elevating leg supports).
- Strengthened suspension for stability.
- Sturdy arm supports and leg supports.
- Wider footplates.
- Corpus seating for the ultimate in comfort.

Size Width 19" - 25" x Depth 19" - 25"
SWL 204Kg



Manual Wheelchair

Aspire Evoke 2 HD

Folding Frame

The Evoke 2 HD wheelchair is designed for a user who wants a lightweight, highly adjustable manual wheelchair.

Swing away and removable arm rests. This facilitates user transfers and also makes for easier transport and storage.

Swing in or out, height adjustable and removable foot rests. This facilitates user transfers and allows optimal positioning for the changing needs.

Strong, nylon upholstery with a comfortable cushion and padded backrest, for increased comfort. A padded overlay is available as an accessory.

Depth adjustable seating with rear angle adjustability, for more accurate positioning, customisability and increased user comfort.

Quick release 24" solid rear tyres. The wheels are easily removed to reduce the weight of the chair for transport and storage.

Size Width 22" - 24" x Depth 16" - 19"
SWL 180Kg



Aspire Evoke 2 HD

Manual Handling

Aspire A320 Aluminium Bariatric Lifter

The Aspire A320 Aluminium Bariatric Lifter will exceed expectations even in the most demanding clinical environments. Its broad base and ample patient clearance compliments the substantial safe working load.

SWL 320Kg



Aidacare Aspire Inline Weigh Scale

Designed to be fitted on Aidacare Aspire Patient Lifters.

Allows for an accurate weight measurement to be taken by a carer when a patient is suspended in a sling.

Easy to read LCD display.

SWL 320Kg



Slings

The Aspire range includes general purpose and bariatric slings with or without head support, hygiene slings for toileting and dressing, all day slings for those who need to remain in their sling for extended periods, and standing slings for patients who can bear weight.

Size S, M, L, XL



ATM-Airpod®

The world's most advanced patient lateral transfer and repositioning system using air flow technology.

Lateral patient transfer and repositioning is a rigorous and potentially injurious process to care-givers and patients. Airpod offers safe patient handling and transfers as well as to help prevent workplace injuries for caregivers.

The ATM-Airpod is engineered to withstand multiple uses under the most difficult transfers for the widest variety of patient profiles.

Our patented technology automatically profiles and self adjusts to the individual patient. This provides a balanced redistribution of the patient's body weight for the safest and easiest patient transfer.

ATM-AirPod features

- HF welding for excellent durability.
- Self regulating and self adjusting.
- Superior hygiene and infection control.
- Advanced pressure redistribution.
- Anti static and MRI translucent.
- Compatible with existing systems.

SWL 400Kg



K Care MAXI Walking Tutor

A rigid, fully welded, four wheeled, mobile walking frame with adjustable height forearm supports with upright foam padded handles. The supports also adjust in cuff length and angle. The frame has swivel front castors and rear wheels or optional glides.


As an option , a padded chest and arm rest is available.

Available in one size.

Size 75cm (between armrest)
SWL 250Kg



TO BE LAUNCHED in 2020

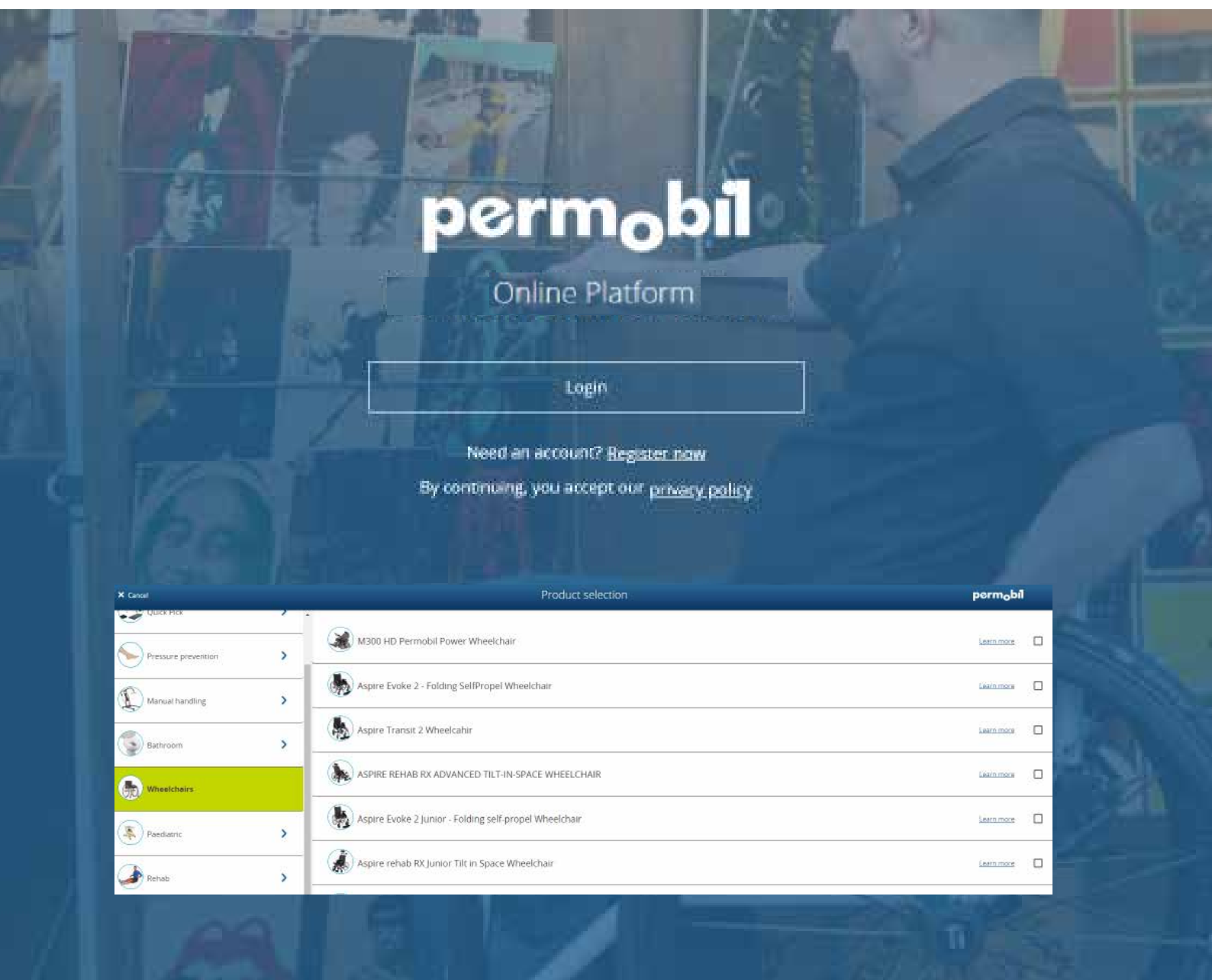

permobil

Online Platform

Login

Need an account? [Register now](#)

By continuing, you accept our [privacy policy](#)

Product selection

		permobil
Quick Pick	M300 HD Permobil Power Wheelchair	Learn more <input type="checkbox"/>
Pressure prevention	Aspire Evoke 2 - Folding SelfPropel Wheelchair	Learn more <input type="checkbox"/>
Manual handling	Aspire Transit 2 Wheelchair	Learn more <input type="checkbox"/>
Bathroom	ASPIRE REHAB RX ADVANCED TILT-IN-SPACE WHEELCHAIR	Learn more <input type="checkbox"/>
Wheelchairs	Aspire Evoke 2 Junior - Folding self-propel Wheelchair	Learn more <input type="checkbox"/>
Paediatric	Aspire rehab RX Junior Tilt in Space Wheelchair	Learn more <input type="checkbox"/>
Rehab		